

(11) **EP 1 118 268 A1**

(12) **EUROPEAN PATENT APPLICATION**

(43) Date of publication:  
 25.07.2001 Bulletin 2001/30

(51) Int Cl.7: **A01N 1/02**

(21) Application number: **00100720.2**

(22) Date of filing: **14.01.2000**

(84) Designated Contracting States:  
**AT BE CH CY DE DK ES FI FR GB GR IE IT LI LU**  
**MC NL PT SE**  
 Designated Extension States:  
**AL LT LV MK RO SI**

(71) Applicant: **ARTEMIS Pharmaceuticals GmbH**  
**51063 Köln (DE)**

(72) Inventors:  
 • **Walderich, Brigitte, Dr.**  
**72076 Tübingen (DE)**

• **Nordin, Renate**  
**72072 Tübingen (DE)**

(74) Representative:  
**Helbing, Jörg, Dr.Dipl.-Chem. et al**  
**Patentanwälte**  
**von Kreisler-Selting-Werner,**  
**Deichmannhaus am Hauptbahnhof**  
**50667 Köln (DE)**

Remarks:

Amended claims in accordance with Rule 86 (2)  
 EPC.

(54) **Method for cryoconservation of zebrafish sperm**

(57) The present invention provides a method for cryoconservation of Zebrafish sperm, and samples of Zebrafish sperm obtainable by said method.

**EP 1 118 268 A1**

## Description

**[0001]** The present invention provides a method for cryoconservation of Zebrafish sperm, and samples of Zebrafish sperm obtainable by said method.

**[0002]** Zebrafish sperm can be prepared by two different methods for cryoconservation and In Vitro Fertilization (IVF) of Zebrafish eggs:

1. Sperm can be squeezed directly from male fish and frozen in cryoprotectant.
2. Sperm can be isolated from prepared fish testis and frozen in cryoprotectant. (Ransom & Zon, Methods in Cell Biology, Vol. 60, p. 365- 372, (1999)). Said method comprises homogenizing the testis by using a Kontes pellet pestle.

Method 1 is established and functions, however it has several disadvantages:

- About 50% of all fish give no sperm sample by using the squeezing procedure.
- Such fish can be squeezed again after 2 weeks: This needs space and time.
- The recovered sperm has a volume of 0,5 - 2  $\mu$ l, resulting in only one aliquot of sperm per squeezing, i.e. the procedure has to be repeated several times.
- The squeezed sperm can differ very much in quality yielding fertilization rates between 0 - 50%.

Method 2 has been described to function in a protocol obtained from Ransom & Zon, 1999, however we have been not able to reproduce it. This method would have several advantages in comparison to the squeezing method:

- Sperm samples can be obtained from every single fish.
- Several sperm aliquots can be frozen from one testis preparation: time and space saving.

**[0003]** Starting from the method of Ransom & Zon, (1999), a reliable method to freeze Zebrafish sperm prepared from testis was established. This method gives at least 2 sperm samples from one testis preparation. The new method can be done very quickly: 5 minutes per fish are needed for the sperm preparation. As demonstrated in 23 independent experiments, fertilization rates between 10 - 50% were observed, using fish between 4 and 12 months old.

**[0004]** The present invention thus provides a method for cryoconservation of Zebrafish sperm comprising

- (a) removing the testis from the body cavity of a male Zebrafish;
- (b) rinsing the testis isolated in step (a) with a cryoprotectant to yield a sperm suspension; and
- (c) freezing said sperm suspension,

wherein the amount of cryoprotectant used in step (b) is such that the volume of the sperm suspension be less than 30  $\mu$ l.

**[0005]** The present invention further provides samples of Zebrafish sperm obtained by the above method.

## Description on the figures

**[0006]** Fig. 1 shows the variation of the cryoprotectant and its concentration on the fertilization rate.

**[0007]** Fig. 2 shows the influence of the aliquot size on the fertilization rate.

**[0008]** Fig. 3 shows the influence of the storage time.

**[0009]** Fig. 4 shows the influence of the age of the male fish on the fertilization.

**[0010]** Fig. 5 shows the influence of the size of the fish on the fertilization rate.

**[0011]** In the following the sperm preparation and cryoconservation is described in more detail.

Male Zebrafish are anaesthetized, dried with a tissue paper, decapitated and the testis is removed from the body cavity. It is important to dry the fish and the equipment used very well, as water may activate the sperm. Activated sperm cannot be used for successful IVFs.

The isolated testis are rinsed very carefully with about 20  $\mu$ l of cryoprotectant having a temperature from 0 to 10°C, preferably 4°C, using a mikrotip.

**[0012]** The suspension obtained should have a volume of at least 20  $\mu$ l but less than 30  $\mu$ l per fish.

The cryoprotectant preferably consists of 1 x Ginzburg Ringer Solution containing 15% wt.% milkpowder and 15% methanol vol.%.

The sperm suspension can be divided into 2 to 5 aliquots, preferably into 5  $\mu$ l portions in glass capillaries.

For freezing the capillaries are cooled down in cryovials and plastic tubes to -70°C for at least 30 minutes and are then transferred in the cryovial into liquid nitrogen, i.e., are cooled down to -196°C (storage temperature).

**[0013]** The *in vitro* Fertilization is performed as follows:

- Female Zebrafish are anaesthetized and dried very carefully.
- The eggs are squeezed by pressing the belly of the female towards the anal fin.
- A capillary with frozen sperm is thawed, mixed with 70  $\mu$ l of Hanks final solution and given to the eggs. Incubation time: 30 seconds at room temperature.
- 750  $\mu$ l fish water is added to the eggs and incubated 2 minutes at room temperature.
- 9 ml of water are added and the eggs are incubated at 28°C for 6 hours.
- The fertilized eggs are transferred into embryo - medium E3.

**[0014]** The following experiments were done: The method to isolate a sperm suspension from testis was altered. Ransom and Zon (1999) describe to homoge-

nize the testis by using a Kontes pellet pestle. Since this method didn't work the sperm were isolated from the testis by flushing the testis very carefully with the cryopreservation solution on ice using different micropipette tips. The most successful method was to use 10 µl cristal tips.

**[0015]** Dimethylsulfoxide (DMSO) was tried instead of methanol, as most cryopreservation methods are using DMSO as a cryoprotectant. Methanol was more efficient compared to DMSO, as shown in the fertilization rates (see Fig. 1).

**[0016]** Different concentrations of DMSO and Methanol were tried: 10%, 15% and 20% of Methanol worked, 15% Methanol gave the best results. 7,5% and 10% of DMSO worked as well, however the fertilization rates obtained were lower than with 15% Methanol (see Fig. 1).

**[0017]** Different volumes of the sperm suspension were tested: 10, 20, 30, 50, 100 and 500 µl. It was shown that sperm samples > 30 µl gave low or negative results in the IVF. The optimal volume of 20 µl can be divided into 4 aliquots a 5 µl.

**[0018]** The cooling down protocol itself was tried to change: Ransom and Zon described that the sperm has to be frozen on crushed CO<sub>2</sub> for 30 minutes and then in liquid nitrogen. This method was shown to be the most efficient compared to freezing the sperm samples directly into liquid nitrogen or cooling down 1°C per minute. These experiments are described in more detail below.

#### Examples

**[0019]** Material used:

- Milk powder (Carnation Nonfat Dry Milk, Nestle)
- Methanol p.a. (Merck, Darmstadt)
- DMSO Hybrimax (Sigma, München)
- RPMI 1640 (Sigma, München)
- Foetal calf serum (Life Technologies, Karlsruhe)
- cristal tips (Eppendorf, Hamburg)
- cryovials, 2 ml (Roth, Karlsruhe)
- 50 ml plastic tubes (screw caps) (Greiner, Nürtingen)
- MESAB Solution (0,4% ethyl-m-aminobenzoate methanesulfonate + 1% Na<sub>2</sub>HPO<sub>4</sub> x 2 H<sub>2</sub>O, pH 7,2)
- 1 x Ginzburg-Ringer-Solution (Ransom & Zon, 1999)
- Cryoprotectant solution: 1 x Ginzburg-Ringer-Solution is prepared and frozen in 4,5 ml aliquots. 15% milk powder (w/v) is resolved and 15% Methanol (v/v) is added.  
The cryoprotectant is cooled down on ice and used freshly prepared.
- Hanks Final Solution (Ransom & Zon, 1999)
- Embryo medium E3 (Haffter et al., 1996 Development 123, p. 1- 36)

#### Example 1

**[0020]** It was tried to reproduce the optimal freezing conditions as described above in 23 independent experiments. In comparison to these tests, the experiments using 10% DMSO instead of 15% Methanol are shown in Fig. 1. One column represents the fertilization rates of at least 6 independent IVFs. The experiments were carried out as described below.

- Male fish (4 - 12 months old) were anaesthized in 4 ml of the MESAB solution per 100 ml fish water, dried very carefully using a tissue paper and decapitated.
- The testis were dissected and transferred into an Eppendorf cup filled with 20 µl of icecold cryoprotectant solution. The cryoprotection solution consisted of 1 x Ginzburg Ringer Solution + 15% Milk powder (w/v) + 15% Methanol (v/v) or of 1 x Ginzburg Ringer Solution + 15% Milk powder (w/v) + 10% DMSO (v/v) and was prepared freshly. The testis were rinsed with the icecold solution using a micropipette with a cristal tip.
- Aliquots of 10, 20 and 30 µl were filled in capillaries and cryotubes and frozen in a 50 ml plastic tube on dry ice for 30 minutes.
- The samples in the cryovials were transferred into liquid nitrogen.
- After one week different sperm samples were used for IVF experiments.
- Female Zebrafish were anaesthized with the MESAB solution (see above) and dried very carefully using a tissue paper.
- The eggs were squeezed by pressing the belly of the fish towards the anal fin into a petri dish.
- 1 aliquot of thawed sperm was resuspended in 70 µl of Hank's final solution, added to the eggs and incubated 30 seconds at room temperature.
- 750 µl of fish water were added and incubated 2 minutes at room temperature.
- 9 ml of fish water were added and the eggs were incubated 6 hours at 28°C.
- Fertilized eggs were separated and counted to calculate the fertilization rate. They were transferred into E3 medium. The unfertile eggs were counted and discarded.

Result: It was shown that under these conditions fertilization rates between 2 and 20% were achieved. In some cases fertilization rates up to 50% and more could be found. The experiments using 15% Methanol as a cryoprotectant gave the highest fertilization rates of 20%. 10% DMSO as cryoprotectant yielded fertilization rates lower than 3,5%. It can also be seen that preparation volumes higher than 20 µl gave lower fertilization rates.

Example 2

**[0021]** 20 µl sperm solution, prepared as described in example 1 with 15% Methanol, were divided into several portions: 1 x 20 µl, 2 x 10 µl, 4 x 5 µl, 10 x 2 µl. These sperm samples were frozen and thawed and used for IVF as described above. 6 independent experiments were done.

Results: 20 µl of sperm solution can be divided into 5 µl portions with a similar fertilization rate as for 10 or 20 µl portions (Fig. 2). 2 µl aliquots show significantly lower fertilization rates.

Example 3

**[0022]** The sperm samples were prepared from 12 month old fish as described above in 10 µl portions and were thawed at different timepoints to show how long they can be stored in liquid nitrogen without losing activity.

Result: We have thawed 4 samples of different fish after 1 week, 2 weeks and 1 month of storage. In these experiments we could not observe lower fertilization rates. (Fig. 3)

Example 4

**[0023]** Sperm samples were prepared as described above from 4 month and 12 month old fish and frozen in 10 µl sperm samples. These were used for IVF in 6 independent experiments as described.

Result: It could be shown that sperm of 12 month old fish yield higher fertilization rates (ca. 42%) than 4 month old fish (ca. 10%). Therefore older fish should be taken for testis preparation rather than younger fish. (Fig. 4)

Example 5

**[0024]** Sperm samples were prepared as described above from 4 month old fish which differed in size. The samples were frozen as 10 µl aliquots and frozen and thawed and used for IVF as described above. 6 samples of big fish were compared to 10 samples of small fish.

Results: The samples of small fish gave lower fertilization rates (ca. 9%) compared to the big fish (ca. 14%). This means that for sperm preparation the fish should be well fed and should be as big as possible to the timepoint where sperm is isolated (Fig. 5)

Claims

1. Method for cryoconservation of Zebrafish sperm comprising

(a) removing the testis from the body cavity of a male Zebrafish;

(b) rinsing the testis isolated in step (a) with a cryoprotectant to yield a sperm suspension; and  
(c) freezing said sperm suspension,

wherein the amount of cryoprotectant used in step (b) is such that the volume of the resulting sperm suspension is less than 30 µl.

2. The method of claim 1, wherein prior to step (a) the male Zebrafish is anaesthetized, dried and/or decapitated.

3. The method of claim 1 or 2, where

- the amount of cryoprotectant used in step (b) is from 10 to 25 µl, preferably 17 to 22 and/or
- the cryoprotectant is a modified 1 x Ginzburg-Ringer-solution containing a polar organic solvent; and/or
- the temperature of the cryoprotectant is from 0 to 10°C, preferably from 2 to 5°C.

4. The method of claim 3, wherein the modified 1 x Ginzburg-Ringer-solution contains to 20 vol. %, preferably 10 to 17 vol.% of the polar organic solvent and may further contain 5 to 25 wt.-%, preferably 10 to 20 wt.-% of a concentrated protein.

5. The method of claim 3 or 4, wherein the polar organic solvent is selected from methanol, DMSO and ethanol, and preferably is methanol.

6. The method of claim 4, wherein the concentrated protein is selected from milk powder, soya and bovine serum albumin and preferably is milk powder.

7. The method according to any one of claims 1 to 6, wherein prior to step (c) the sperm suspension is divided into 2 to 5 aliquots, preferably into 5 µl portions.

8. The method of claim 7, wherein the portions are kept in glass capillaries.

9. The method according to any one of claims 1 to 8, where in step (c) the sperm suspension is kept at -70°C (± 10°C) for at least 30 min and is then cooled down to the temperature of liquid nitrogen.

10. Samples of cryoconserved Zebrafish sperms obtainable by the method of claims 1 to 9.

Amended claims in accordance with Rule 86(2) EPC.

1. Method for cryoconservation of Zebrafish sperm comprising

(a) removing the testis from the body cavity of a male Zebrafish;

(b) rinsing the testis isolated in step (a) with a cryoprotectant to yield a sperm suspension; and

(c) freezing said sperm suspension,

wherein the amount of cryoprotectant used in step (b) is such that the volume of the resulting sperm suspension is less than 30  $\mu$ l.

2. The method of claim 1, wherein prior to step (a) the male Zebrafish is anaesthetized, dried and/or decapitated.

3. The method of claim 1 or 2, where

- the amount of cryoprotectant used in step (b) is from 10 to 25  $\mu$ l, preferably 17 to 22  $\mu$ l; and/or
- the cryoprotectant is a modified 1 x Ginzburg-Ringer-solution containing a polar organic solvent; and/or
- the temperature of the cryoprotectant is from 0 to 10°C, preferably from 2 to 5°C.

4. The method of claim 3, wherein the modified 1 x Ginzburg-Ringer-solution contains 5 to 20 vol. %, preferably 10 to 17 vol.% of the polar organic solvent and may further contain 5 to 25 wt.-%, preferably 10 to 20 wt.-% of a concentrated protein.

5. The method of claim 3 or 4, wherein the polar organic solvent is selected from methanol, DMSO and ethanol, and preferably is methanol.


6. The method of claim 4, wherein the concentrated protein is selected from milk powder, soya and bovine serum albumin and preferably is milk powder.


7. The method according to any one of claims 1 to 6, wherein prior to step (c) the sperm suspension is divided into 2 to 5 aliquots, preferably into 5  $\mu$ l portions.


8. The method of claim 7, wherein the portions are kept in glass capillaries.

9. The method according to any one of claims 1 to 8, where in step (c) the sperm suspension is kept at -70°C ( $\pm$  10°C) for at least 30 min and is then cooled down to the temperature of liquid nitrogen.

**Fig. 1: Variation of the cryoprotectant and its concentrations**

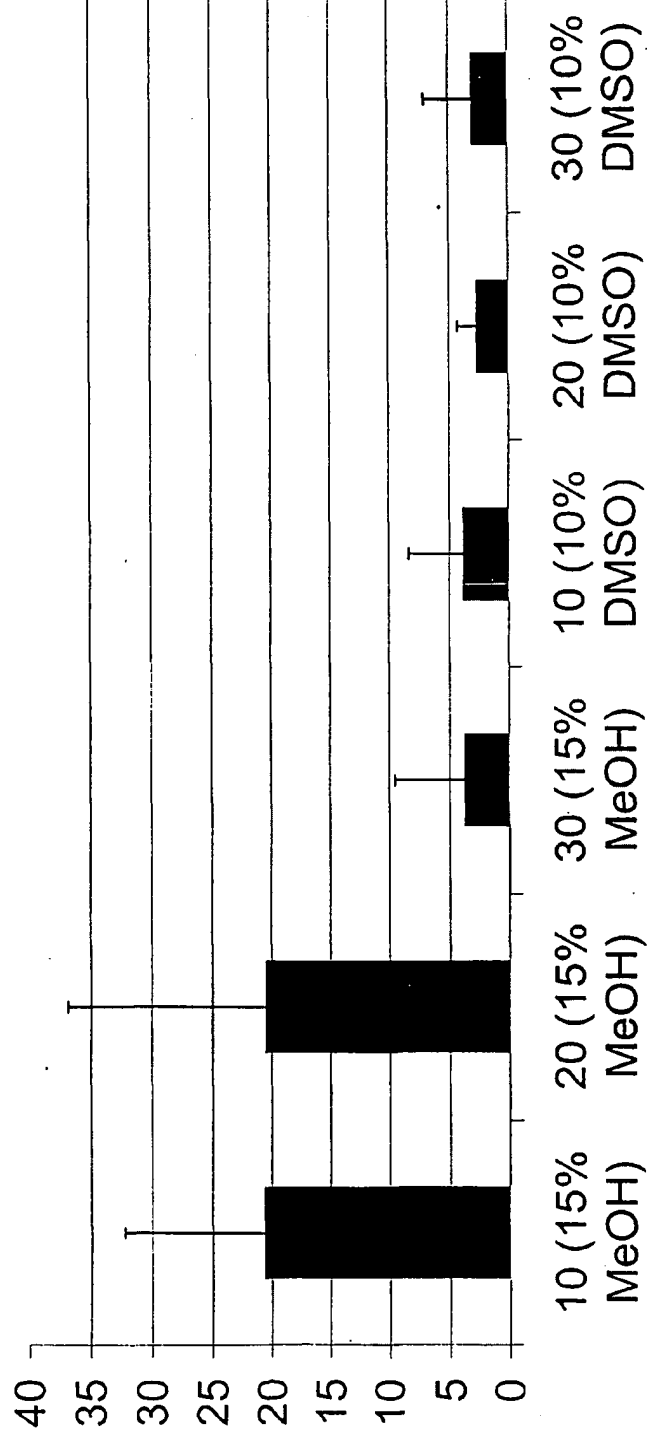


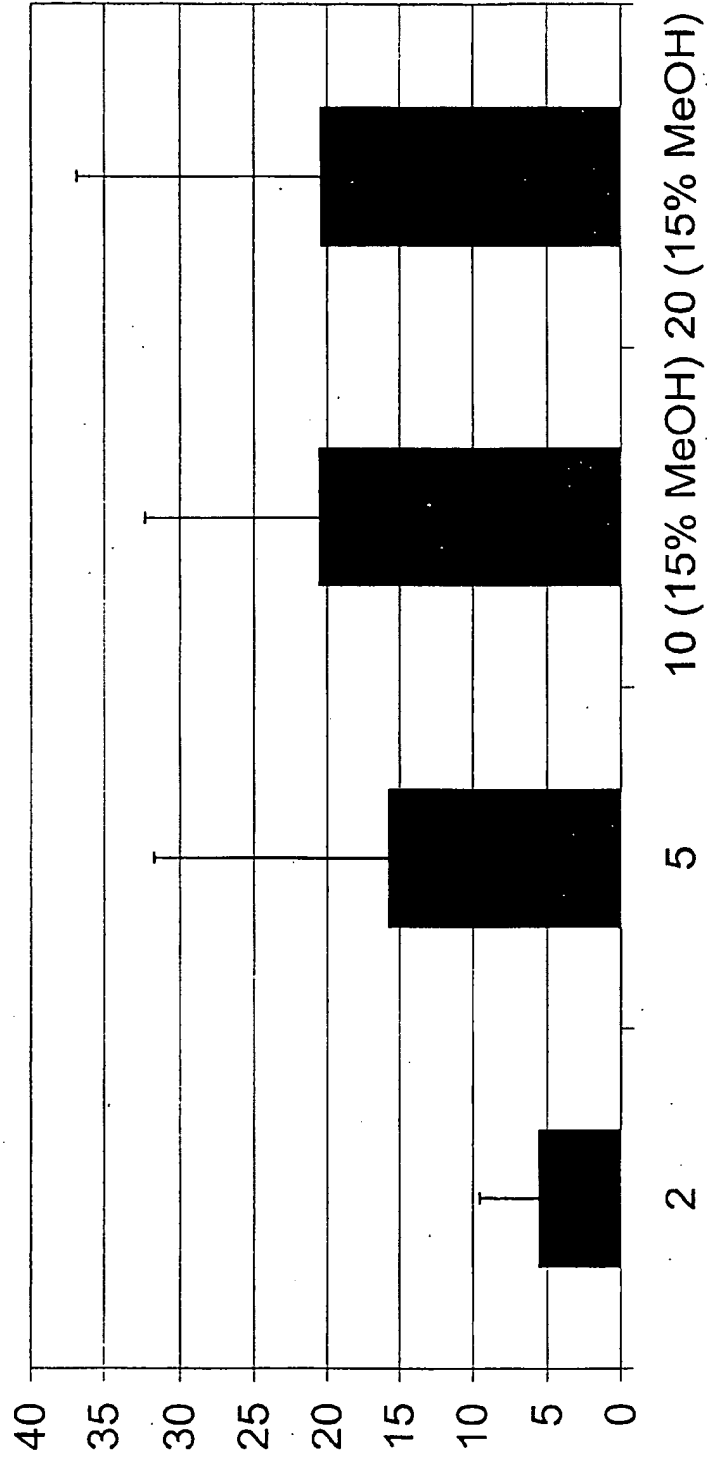
**Fig. 2: Aliquots**

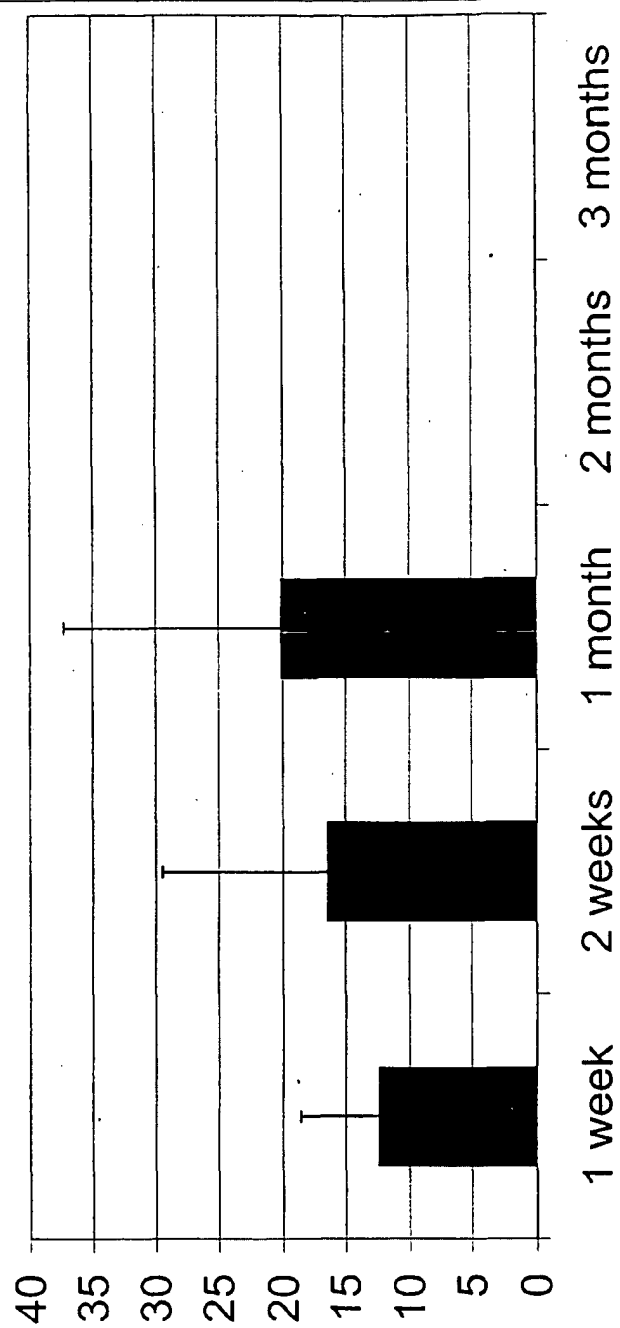**Fig. 3: Long term storage**



**Fig. 4: 4 month/12 month old fish**

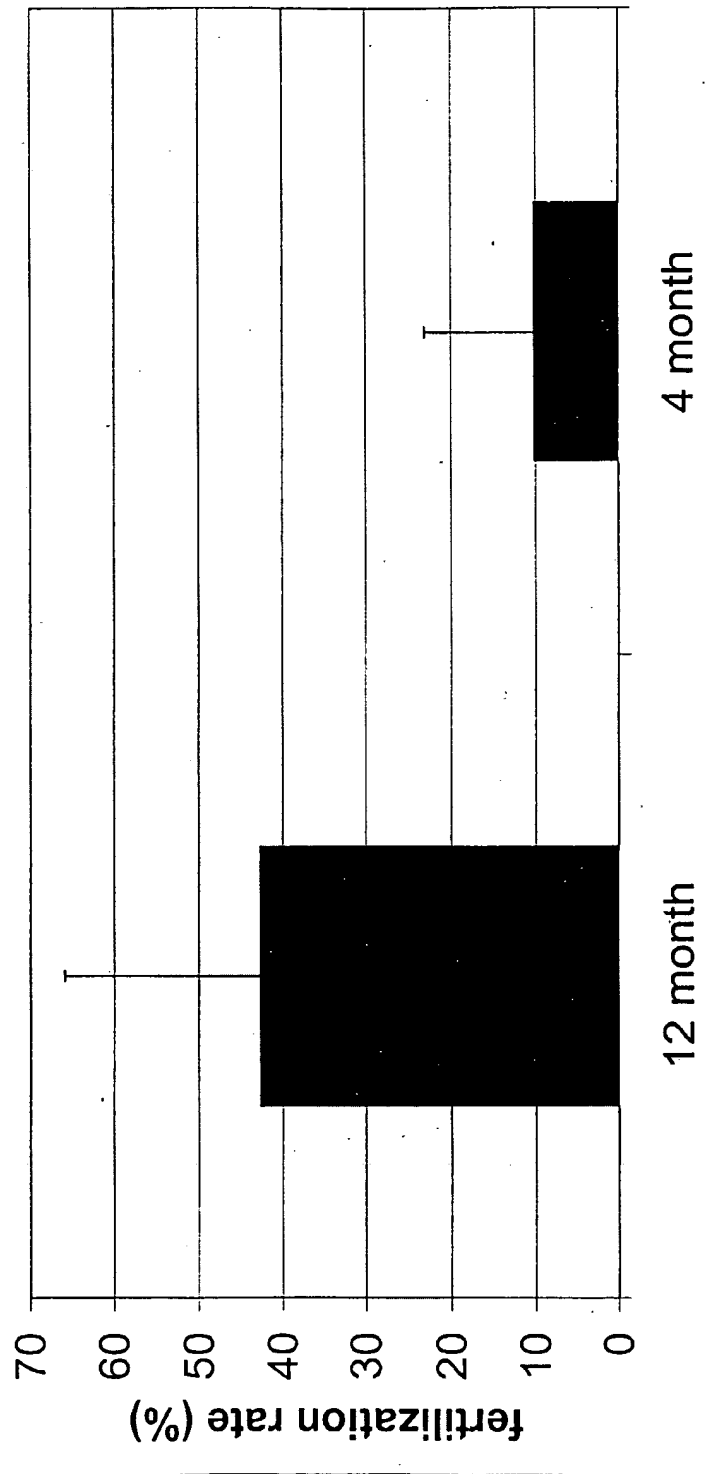

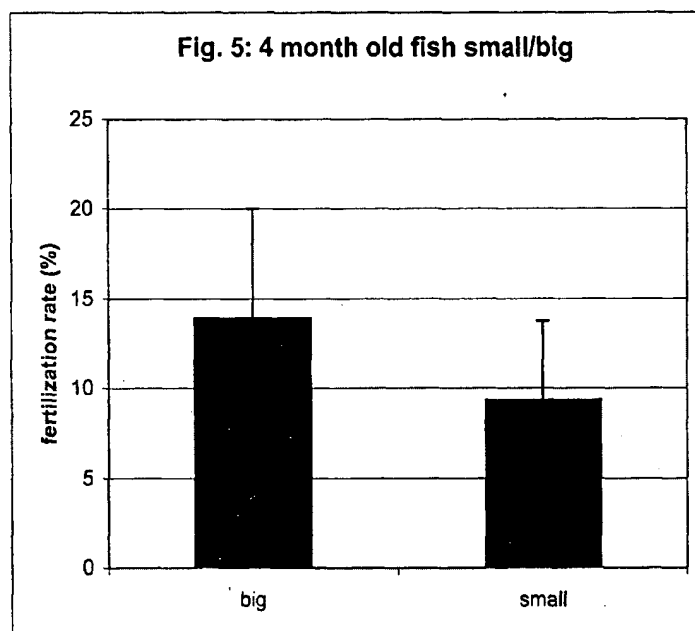



European Patent  
Office

# EUROPEAN SEARCH REPORT

Application Number  
EP 00 10 0720

DOCUMENTS CONSIDERED TO BE RELEVANT			
Category	Citation of document with indication, where appropriate, of relevant passages	Relevant to claim	CLASSIFICATION OF THE APPLICATION (Int.Cl.7)
D,A	D. G. RANSOM & L. I. ZON: "Collection, storage and use of Zebrafish sperm" METHODS IN CELL BIOLOGY, vol. 60, 1999, pages 365-372, XP000909881	1-9	A01N1/02
X	* page 369, paragraph B *	10	
A	B. HARVEY, R. N. KELLEY & M. J. ASHWOOD-SMITH: "Cryopreservation of zebra fish spermatozoa using methanol" CAN. J. ZOOLOGY, vol. 60, no. 8, 1982, pages 1867-1870, XP000909956	1-9	
X	* page 1867, column 2, paragraphs 2,3 *	10	
A	CH. WALKER: "Freezing sperm" ZEBRAFISH SCIENCE MONITOR, 'Online! vol. 1, no. 3, 1991, pages 4-5, XP002139737	1-9	
	Retrieved from the Internet: <URL:http://zfin.org/zf_info/monitor/mon.html> 'retrieved on 2000-05-30!		TECHNICAL FIELDS SEARCHED (Int.Cl.7)
X	* the whole document *	10	A01N
A	& M. WESTERFIELD: "The Zebrafish Book, 3rd. ed." 'Online! 1995, UNIV. OREGON PRESS, EUGENE, OREGON	1-9	
	Retrieved from the Internet: <URL: http://zfin.org/zf_info/zfbook/zfbk.html> 'retrieved on 2000-05-30!		
X	See: chapter 2.8: Breeding. Embryo production by in vitro fertilization.	10	
The present search report has been drawn up for all claims			
Place of search MUNICH		Date of completion of the search 8 June 2000	Examiner: Klaver, J
<p>CATEGORY OF CITED DOCUMENTS</p> <p>X : particularly relevant if taken alone Y : particularly relevant if combined with another document of the same category A : technological background O : non-written disclosure P : intermediate document</p> <p>T : theory or principle underlying the invention E : earlier patent document, but published on, or after the filing date D : document cited in the application L : document cited for other reasons &amp; : member of the same patent family, corresponding document</p>			

EPO FORM 1502 01.02 (P0401)